


QC Report – Sample GSM8461135

Date:2026-05-21; Library: ewastools

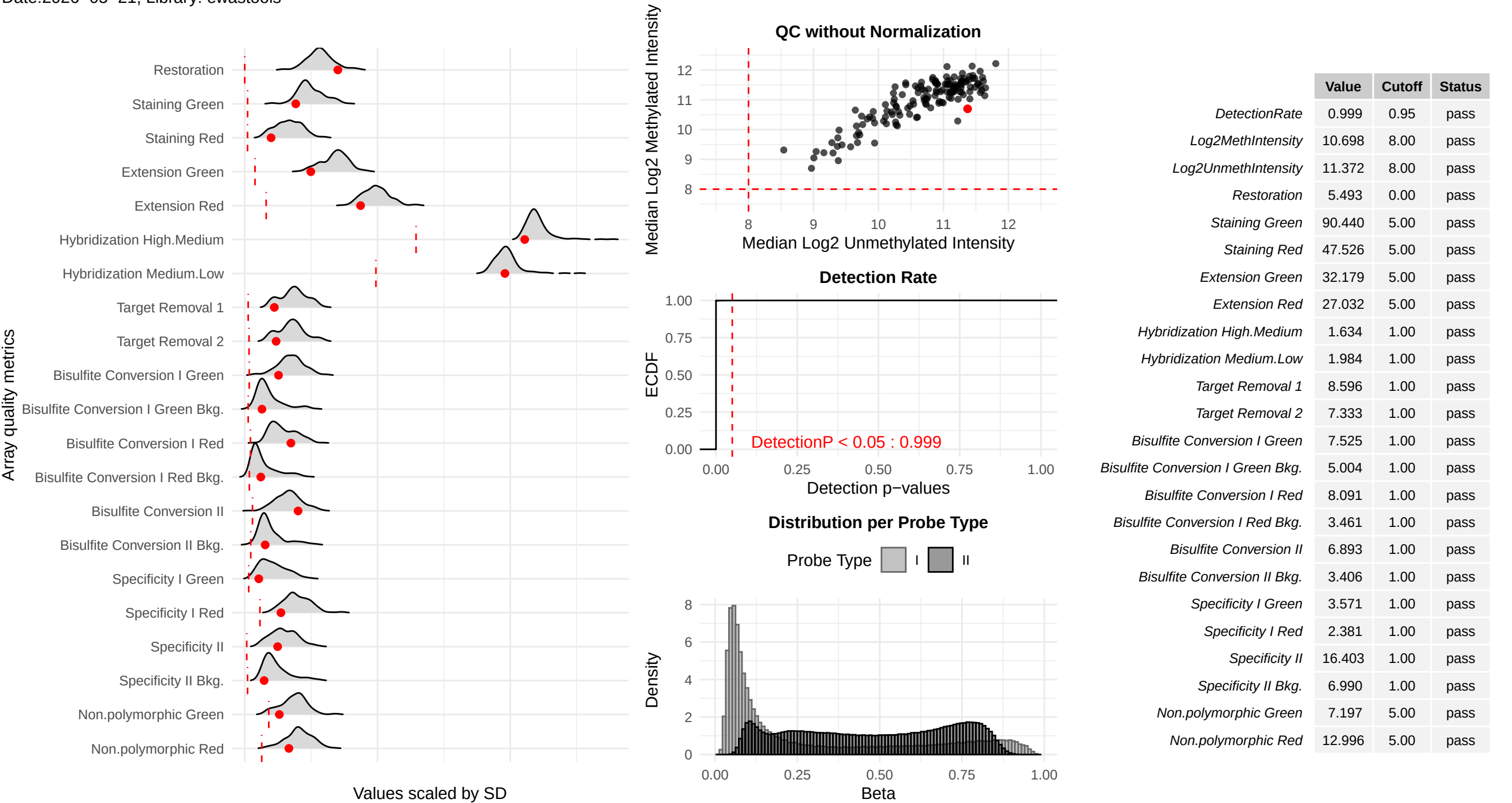

On the left, the figure shows 21 quality metrics derived from the control probes on the Illumina Methylation Arrays. Red dots indicate the values for the analyzed sample, while the surrounding distributions reveal the variability observed in our database. Red vertical lines mark the recommended thresholds. Overall array quality is assessed using three primary criteria:

- (1) the proportion of CpG probes with a detection p-value below 0.05,
- (2) the log2 median intensity of methylated signal,
- (3) the log2 median intensity of unmethylated signal.

A sample is considered acceptable for methylation classification when all three of these benchmarks meet or exceed their thresholds. Although slight deviations in other quality metrics may not impact classification accuracy, they might point to issues during laboratory protocols, hence, simultaneous failures in multiple metrics should be investigated further.